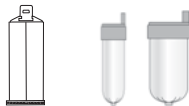
# Monopren® transfer

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35713 Eschenburg, Germany

[www.kettenbach.com](http://www.kettenbach.com)

**Vinyl Polysiloxane Impression Material**

**ISO 4823**

Made in Germany  
37203/1311

For professional use only.

# Vinyl Polysiloxane Impression Material ISO 4823

## Product Description

Monopren® transfer, the medium-viscosity monophase impression material, offers significant advantages for taking impressions during combined prosthetic work and in implantology. Thanks to the initial hydrophilicity, optimum hydrophilic properties appear after just 0.5 seconds. These properties promote immediate inflow to the primary crowns, cast support or preparation limits even in the moist oral environment. Other important criteria include good setting characteristics, the appropriate Shore A hardness and thixotropy. As a monophase material, it is an all-purpose solution.

## Indications/Techniques

- Medium viscosity, monophase impression material for taking impressions over fixed/removable restorations and implants (transferring impression posts and bridge components).
- Medium viscosity impression material for functional impressions.

- Medium viscosity, monophase impression material for use when fabricating crown and bridgework or inlays.
- Medium viscosity impression material for reline impressions.
- Syringeable impression material for use in the simultaneous mixing technique as well as the putty-wash and triple tray techniques.
- Medium viscosity, monophase impression material for transferring root posts when fabricating posts and cores indirectly.
- Copper ring impressions.

## Warnings

Do not use Monopren® transfer impression materials as a temporary reliner.

Do not use Monopren® transfer impression materials with condensation-curing silicones, polyether, Vinylsiloxanether® or polysulfide materials.

## Cautions

Do not use after expiration date.

Do not leave any residual impression material in the sulcus or oral cavity.

Do not swallow impression material! If swallowed: consult a medical doctor if problems arise or persist.

Avoid contact with eyes. If accidental contact with the eyes occurs, rinse immediately and thoroughly with an eye wash or water. Consult a medical doctor if problems arise or persist.

The product does not normally cause allergic reactions. However, for sensitive persons, an allergic reaction to the product cannot be ruled out.

Use of products containing active sulfur, aluminium chloride or nitrogen compounds (retraction cords

containing ferric sulfate, polysulfide impression materials, etc.) in conjunction with this product will interfere with the setting reaction of the vinyl polysiloxane material. Use of the materials requires the area to be rinsed thoroughly to remove all residue, before an impression is taken. Do not use latex gloves.

Before fitting the dynamic mixer, dispense material until equal amounts of base and catalyst appear; wipe excess. Firmly place dynamic mixer and lower locking lever.

When taking impressions of areas with severe undercut and wide interdental spaces, use standard block-out techniques.

When using a custom impression tray, ensure that sufficient space remains between the side of the tray and the teeth/jaw.

Do not allow the material to enter the sewer or water system to avoid environmental contamination.

Avoid contact with clothing, since the material cannot be removed by dry cleaning.

## Product Overview

Product name	ISO 4823	Consistency (approx.) mm	Mixing ratio and total content volume per unit	Mixing technique	Mixing element
Monopren® transfer	Type 2, Medium-bodied	35	5:1, 362 ml in foil bags	Plug & Press® Dispenser, Sympress or other automatic dispensing and mixing unit	dynamic mixer
Monopren® transfer	Type 2, Medium-bodied	35	1:1, 50 ml cartridge	Applyfix® 4 dispensing gun DS-50 1:1/2:1	static green mixing tip MB Ø 6.5 mm

## Technical Data

Product name	Mixing ratio	Working time at 23° C/ 74° F ≤	Working time at 35° C/ 95° F ≤	Intraoral setting time at 35° C/ 95° F ≥	Total setting time* ≥	Hardness (approx.) Shore	Linear dimensional change (approx.) %	Elastic-Recovery test (approx.) %	Strain in compression (approx.) %
Monopren® transfer	5:1	2 minutes	1 minute	2 minutes	4 minutes	A 60	- 0.20	99.7	3.5
Monopren® transfer	1:1	2 minutes	1 minute	2 minutes	4 minutes	A 60	- 0.20	99.7	3.5

\*Total setting time (removal time from mouth) from start of mix.

## Note

To ensure optimal impressions, the temperature of the material should not deviate from 23° C (74° F) before applying. Otherwise, working and setting times will be affected.

## Impression tray: Preparation and adhesives

In principle, all common impression trays can be used if a relevant dynamic pressure is guaranteed. When retention is not sufficient, ensure a strong bond to the impression material, by brushing the impression tray with a thin film of Panasil® adhesive for addition-curing silicone prior to loading the tray with impression material. Allow to dry per manufacturer's instructions.

## Preparation of the Plug & Press® foil bag system

### Foil bag material for use in most automatic dispensing and mixing systems

Before using the base material for the first time, remove the white safety pin from the activating head

of the large foil bag by turning it in the direction of the arrow and pulling it out (Figure 1, Page 80/81).

Place the large foil bag with the integrated activating head in the cartridge body. Ensure that the notches on the activating head and cartridge case are aligned (Figure 2, Page 80/81).

Firmly press the activating head into the final position on the cartridge body. The foil bag is automatically pierced by a pin when the activating head is pressed down (Figure 3, Page 80/81).

Follow the same procedure as in Figure 2 and 3 for the catalyst material. (Note: There is no safety pin on the catalyst foil bag.)

It is recommended to use only Kettenbach dynamic mixers to ensure optimum performance.

## Application using the Plug & Press® Dispenser in combination with the Kettenbach's dynamic mixer

- Move the plungers to the top by turning the control knob (Figure 4, Page 80/81).
- Insert the cartridge body into the unit (Figure 5, Page 80/81).

- Turn the control knob to move the plungers into the cartridge body, and continue turning until the plungers come into contact with the foil bags (Figure 6, Page 80/81).
- Once there is contact, press one of the two start keys (at the front on top of the unit) to dispense the material (Figure 7, Page 80/81).
- Before fitting the dynamic mixer, dispense material until equal amounts of base and catalyst appear (Figure 7, Page 80/81); wipe excess. Firmly place dynamic mixer (Figure 8, Page 80/81) and lower locking lever.

Load the tray with the required amount of material. Hold the impression tray at an angle and press lightly against the tray. Leave the dynamic mixer in the material while dispensing (Figure 9, Page 80/81). Leave the filled dynamic mixer on the cartridge body as a seal.

Before next use, release the locking clip to remove the used dynamic mixer and check that the outlets in the activating heads are not blocked. Fit a new dynamic mixer, lower locking lever and continue as usual.

When finished, simply remove the empty foil bags and activating heads from the cartridge body and discard them. The cartridge body is reusable (Figure 10, Page 80/81). The activating heads are disposable.

If not using a Kettenbach Plug & Press® Dispenser, please follow the manufacturer's instructions for your automatic dispensing and mixing system.

## **Application using the Sympress Dispenser in combination with the Kettenbach's dynamic mixers**

- Allow the plungers (with cover closed) to return to the original position by pressing the plunger return button (Figure 1, Page 82/83). Next, open the cover (Figure 2, Page 82/83).
- Insert the cartridge body with the foil bags into the device (if foil bag material is already in use, the filled mixer will be on the cartridge body as the closure) (Figures 2 and 3, Page 82/83)
- Mount the dynamic mixer tip in place (if you are already using foil bag material, first remove the mixer, which is being used as the closure and

replace it with the new dynamic mixer tip). (Figure 4, Page 82/83)

- Close the locking lever (Figure 5, Page 82/83)
- Close the device cover (Figure 6, Page 82/83)
- Dispense the material, dispose of the initial approx. 3 cm (until a homogenous mix of base and catalyst appear; see coloring of the final impressions material) (Figures 7 and 8, Page 82/83)
- Load the tray (or application syringe) (Figure 9, Page 82/83). As you do, hold the impression tray at an angle and press lightly against it. Leave the mixer in the material while dispensing to prevent air from getting trapped.
- Leave the filled dynamic mixer in place as a seal.

Before next use, release the locking clip to remove the used dynamic mixing tip and check that the outlets in the activating heads are not blocked. Fit a new dynamic mixer, lower locking lever and continue as usual.

When finished, simply remove the empty foil bags and activating heads from the cartridge body and discard them. The cartridge body is reusable (Figure 10, Page 82/83). The activating heads are disposable.

If you are using an automated dosage and mixing system other than the Sympress Dispenser, please follow the instructions for your respective system.

### **Instructions for use:**

#### **Cartridge material for use in Applyfix® 4 dispensing gun DS-50 1:1/2:1**

Insert the cartridge into the Applyfix® 4 dispensing gun DS-50 1:1/2:1. Ensure that the notches on the cartridge base are pointing down. The clasp will not close if the cartridge is not properly inserted into the dispensing gun.

Remove the cartridge cap. The cap can be replaced after initial use.

Insert plunger into cartridge and dispense a small amount of impression material until equal amounts are dispensed at the same rate.

Install a mixing tip on the cartridge, and turn cap 90 degrees clockwise to lock in place.

Inject the required amount of material, directly into the tray or the preliminary impression, or into a delivery syringe. To apply the material around the prepared teeth, use the delivery syringe or the cartridge with an intraoral tip on the mixer. Check to ensure materials are set prior to mouth removal.

Leave used mixing tip on the cartridge after use or replace the cartridge cap. Before using the cartridge again, remove cartridge cap or remove and discard the used mixing tip.

Check the cartridge openings to ensure no polymerized material is present. Should this be the case, unblock the openings and dispense a small amount of impression material until equal amounts are dispensed at the same rate.

Install a new mixing tip and continue as described above.

The Applyfix® 4 dispensing gun DS-50 1:1/2:1 can be sterilized in the autoclave.

## **Storage conditions:**

Store in a dry place at room temperature; do not expose to sunlight and avoid exposure to extreme temperature fluctuations.

## **Disinfection:**

Impressions can be disinfected using for example, a 2% acidic glutaraldehyde solution. Use disinfection solutions specific for impression materials (e.g., Silosep®). See manufacturer instructions.

## **Electroplating:**

Impressions may be silver-plated with an alkaline-plating solution only.

## **Model preparation:**

For model fabrication, the impression can be cast immediately after disinfection. Impressions can also be poured within two weeks with standard dental stone class IV (e.g., Tewestone® or Tewaterock®).

## **Trademarks**

Kettenbach®, Applyfix® 4, Monopren® transfer, Panasil®, Silosept®, Tewaterock®, Tewestone® and

Plug & Press® Dispenser are registered trademarks of Kettenbach GmbH & Co. KG.

## Note

Selected Kettenbach materials are available in certain markets only.

## Limitation of liability

Except where prohibited by law, Kettenbach GmbH will not be liable for any loss or damage arising from this product, whether direct, indirect, special, incidental or consequential, regardless of the theory asserted, including warranty, contract, negligence, or strict liability.

The information provided for Kettenbach products is based on comprehensive research and experience in application technology.

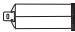
Results are furnished to the best of our knowledge, subject to technical changes within the framework of product development.

However, users must comply with and consider all recommendations and information in connection with any use.

Updated: 28th March 2011

## Legend

Conformity with MDD 93/42/EEC		Use by	
Temperature limitation		Keep dry	
Keep away from sunlight		Catalog number	REF
Batch code		Caution, consult accompanying documents	
Dual system (Germany only)		Diameter	
Cartridge		Static mixing tip MB (single use)	
Mixing tip Type B	MB	Intraoral tip (single use)	
Adhesive		Single use	
Cartridge body		Foil bags	
Dynamic mixer (single use)		Rx only For professional use only.	
Millimeter	mm	Milliliter	ml
Less than or equal to	≤	Greater than or equal to	≥

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